



PRODUCT SPECIFICATION SHEET

Antibiotic Assay Medium A (pH 7.9) (DM023B)

Intended Use

Antibiotic Assay Medium A (pH 7.9) (DM023B) is recommended for the microbiological assay of antibiotics in compliance with BP.

Product Summary and Explanation

Antibiotic Assay media are used in the performance of antibiotic assays. Grove and Randall have elucidated those antibiotic assays and media in their comprehensive treatise on antibiotic assays.⁽¹⁾ The activity (potency) of an antibiotic can be demonstrated under suitable conditions by its inhibitory effect on microorganisms.⁽²⁾ Reduction in antimicrobial activity may reveal changes not demonstrated by chemical methods.⁽²⁾ Antibiotic assays are performed by the cylinder plate method and the turbidimetric "tube" assay. The cylinder plate method, first described by Abraham et al.⁽³⁾ for the assay of penicillin, was later modified by Foster and Woodruff⁽⁴⁾ and by Schmidt and Moyer.⁽⁵⁾ Antibiotic Assay Medium A (pH 7.9) is recommended by BP⁽⁶⁾ and FDA.⁽⁷⁾

Principles of the Procedure

Antibiotic Assay Medium A (pH 7.9) contains a combination of peptone, pancreatic digest of casein, yeast extract and beef extract which provides essential nutrients for the growth of test organisms. Glucose monohydrate provides the carbon source, enhances the growth of test organism. Agar provides excellent medium for antibiotic diffusion and gives well-defined zones of inhibition. Higher pH provides the optimal conditions for activity of antibiotic and also supports the growth of the test organisms.

Cylinder Plate Assay

This method is based on the diffusion of an antibiotic solution from a cylinder placed on the surface of an inoculated agar medium. After incubation, the diameter of a zone of inhibition is observed which depends, in part, on the concentration or activity of the antibiotic. The results depend on critical rates of diffusion of the antibiotic, critical growth rates of the standard organisms and critical minimal inhibitory coefficient levels of each organism. This method is used in the assay of commercial preparations of antibiotics, as well as in the quantitative determination of antibiotics in body fluids, animal feeds and other materials.

Pre-diffusion of antibiotics for 10-20 minutes in the agar by incubating at temperature below the optimal growth temperature for microorganism would facilitate better diffusion of antibiotics followed by incubation of plates for microbial growth.

Note :

Recommended for the microbiological assay of Gentamicin sulphate, Kanamycin monosulphate, Kanamycin acid sulphate, Neomycin sulphate, Netilmycin sulphate, Spiramycin, Streptomycin sulphate, Tylosin, Tylosin tartarate, Vancomycin hydrochloride.

Formula / Liter

Ingredients	Gms / Liter
Peptone	6.00
Pancreatic digest of casein	4.00
Yeast extract	3.00
Beef extract	1.50
Glucose monohydrate	1.00
Agar	15.00





PRODUCT SPECIFICATION SHEET

Final pH: 7.9 ± 0.1 at 25°C
Formula may be adjusted and/or supplemented as required to meet performance specifications

Precautions

1. For Laboratory Use only.
2. IRRITANT. Irritating to eyes, respiratory system, and skin.
3. Freshly prepared plates should be used for antibiotic assays.
4. All conditions in the microbiological assay must be controlled carefully.
5. The use of standard culture medium in the test is one of the important steps for obtaining good results.

Directions

1. Suspend 30.4 grams of the medium (the equivalent weight of dehydrated medium per litre) in one liter of distilled water.
2. Heat to boiling to dissolve the medium completely.
3. Autoclave at 121°C , 15 psi pressure, for 15 minutes / validated cycle.

Quality Control Specifications

Dehydrated Appearance	Cream to yellow coloured homogeneous free flowing powder
Prepared Medium	Light yellow coloured clear to slightly opalescent gel forms in Petri plates
Reaction of % solution	7.9 ± 0.1
Gel Strength	Firm, comparable with 1.5% Agar gel

Expected Cultural Response: Cultural characteristics observed after an incubation at specified temperature for 18-24 hours.

Sr. No.	Organisms	Results to be achieved			
		Inoculum (CFU)	Growth	Recovery	Antibiotics Assayed
1.	<i>Micrococcus luteus</i> ATCC 9341	50-100	good-luxuriant	$\geq 70\%$	Tylosin, Tylosin tartarate (adjust the pH to 8.0 ± 0.1) $32-35^{\circ}\text{C}$
2.	<i>Staphylococcus aureus</i> ATCC 6538p	50-100	good-luxuriant	$\geq 70\%$	Kanamycin monosulphate $30-37^{\circ}\text{C}$ Kanamycin acid sulphate $35-39^{\circ}\text{C}$ Netilmycin sulphate $32-35^{\circ}\text{C}$
3.	<i>Staphylococcus epidermidis</i> ATCC 12228	50-100	good-luxuriant	$\geq 70\%$	Gentamycin sulphate $35-39^{\circ}\text{C}$
4.	<i>Bacillus pumilis</i> NCTC 8241	50-100	good-luxuriant	$\geq 70\%$	Gentamycin sulphate $35-39^{\circ}\text{C}$ Neomycin sulphate $30-37^{\circ}\text{C}$





PRODUCT SPECIFICATION SHEET

5.	<i>Bacillus subtilis</i> ATCC 6633	50-100	good-luxuriant	≥70%	Kanamycin monosulphate 30-37°C, Kanamycin acid sulphate 35-39°C, Neomycin sulphate 30-37°C, Spiramycin 30-32°C, Streptomycin sulphate 30-37°C Vancomycin hydrochloride (adjust the pH to 8.0±0.1) 37-39°C
6.	<i>Bacillus subtilis</i> NCTC 8236	50-100	good-luxuriant	≥70%	Streptomycin sulphate 30-37°C

The organisms listed are the minimum that should be used for quality control testing.

Test Procedure

Preparation of Stock cultures

1. Maintain stock cultures on agar slants and make transfers at 1- or 2-week intervals.
2. Using sterile purified water, saline or Antibiotic Medium No. 3, prepare the inoculum for assay by washing growth from a fresh 24-48 hour agar slant and further dilute the culture to obtain the desired organism concentration.
3. For using as a test organism *Bacillus subtilis*, inoculate the organism on Antibiotic Medium No. 1 and incubate at 37°C for 1 week. Wash spores from the agar surface, and heat the spores at 56°C for 30 minutes. Using sterile purified water wash the spores three times, heat again at 65°C for 30 minutes, and then dilute to the optimal concentration. Inoculum prepared following this method should produce a sharp zone in the assay.
4. For preparing spore suspension of *B. subtilis*, Antibiotic Medium No. 1 modified by the addition of 300 mg manganese sulfate (MnSO₄·H₂O) per liter is used which aids in the sporulation of *B. subtilis*.

Cylinder Plate Assay

1. Use 20 × 100 mm glass or plastic Petri dishes with sufficient depth so that cylinders used in the assay will not be pushed into the medium by the cover.
2. Use stainless steel or porcelain assay cylinders having the following dimensions (± 0.1 mm): 8 mm outside diameter, 6 mm inside diameter and 10 mm long. Clean the cylinders carefully to remove all residues, using an occasional acid bath (i.e., with approximately 2N nitric acid or with chromic acid).
3. Four or six cylinders are generally used per plate, evenly spaced on a 2.8 cm radius.
4. For assuring accurate assays, use a level surface for working to obtain uniformly thick base and seed layers in the Petri dish.





PRODUCT SPECIFICATION SHEET

5. Allow the base layer to solidify and then overlay the seed layer containing a proper concentration of the test organism. The amount of medium in the layers varies for different antibiotics, with most assays specifying a 21 mL base layer and a 4 mL seed layer.
6. In any case, dishes with flat bottoms are required to assure complete coverage of the bottom of the dish when small amounts of base medium are used. Tilt the plate to obtain even coverage of the base layer by the seed layer and allow it to solidify in a level position. Plates should be used the same day as prepared.

Results

1. After incubation the concentration of the antibiotic being assayed is determined by measuring the zone of inhibition obtained, with that of reference standard antibiotic.
2. Refer to appropriate references and specific test procedures.

Storage

Store the sealed bottle containing the dehydrated medium at 10 - 30°C. Once opened and recapped, place container in a low humidity environment at the same storage temperature. Protect from moisture and light.

Expiration

Refer to the expiration date stamped on the container. The dehydrated medium should be discarded if not free flowing, or if the appearance has changed from the original color. Expiry applies to medium in its intact container when stored as directed.

Limitations of the Procedure

1. For identification, organisms must be in pure culture. Morphological, biochemical and/or serological tests should be performed for final identification.
2. Consult appropriate texts for detailed information and recommended procedures.

Packaging

Product Name : Antibiotic Assay Medium A (pH 7.9)

Product Code : DM023B

Available Pack sizes : 100gm / 500gm

References

1. Grove and Randall. 1955. Assay methods of antibiotics. Medical Encyclopedia, Inc. New York, N.Y.
2. United States Pharmacopeial Convention, Inc. 2008. The United States pharmacopeia 31/The national formulary 26, Supp. 1, 8-1-08, online. United States Pharmacopeial Convention, Inc., Rockville, Md.
3. Abraham, Chain, Fletcher, Florey, Gardner, Heatley and Jennings. 1941. *Lancett* *ii*: 177.
4. Foster and Woodruff. 1943. *J. Bacteriol.* *46*:187.
5. Schmidt and Moyer. 1944. *J. Bacteriol.* *47*:199.
6. British Pharmacopoeia, 2011 The Stationery office British Pharmacopoeia.
7. Tests and Methods of Assay of Antibiotics and Antibiotic containing Drugs, FDA, CFR, 1983. Title 21, part 436, Subpart D, Washington, D.C. U.S Government printing office, paragraphs 436, 100-436, 106 pg 242-259 (April 1).

Further Information

For further information please contact your local MICROMASTER Representative.





PRODUCT SPECIFICATION SHEET



MICROMASTER LABORATORIES PRIVATE LIMITED


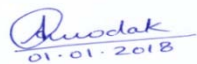

DM023BPSS,QAD/FR/024,Rev.00/01.01.2018

Unit 38/39, Kalpataru Industrial Estate,
Near Runwal Estate, Behind 'R-Mall' ,Ghodbunder Raod,
Thane (W) - 400607. M.S. INDIA.

Ph: +91-22-25895505, 4760, Cell: 9320126789.

Email: micromaster@micromasterlab.com

sales@micromasterlab.com

	Checked By	Approved By
 01.01.2018	 01.01.2018	 01.01.2018
Microbiologist	Head Quality Control	Head Quality Assurance

Disclaimer :

All Products conform exclusively to the information contained in this and other related Micromaster Publications. Users must ensure that the product(s) is appropriate for their application, prior to use. The information published in this publication is based on research and development work carried out in our laboratory and is to the best of our knowledge true and accurate. Micromaster Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are intended for laboratory, diagnostic, research or further manufacturing use only and not for human or animal or therapeutic use, unless otherwise specified. Statements included herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.

