



PRODUCT SPECIFICATION SHEET

Pseudomonas Agar Medium For Detection of Pyocyanin (DM219I)

Intended Use

Pseudomonas Agar Medium For Detection of Pyocyanin (DM219I) is recommended for detection of pyocyanin production by *Pseudomonas* species in compliance with IP.

Product Summary and Explanation

Pseudomonas aeruginosa is one of the most commonly isolated pathogens, and is widely distributed in soil, water and foods. It is the most frequently isolated non-fermentative bacillus in clinical specimens, infusion fluids, disinfectants and cosmetics. The organism causes disease in humans; e.g., ocular infections, burn wound infections and respiratory tract infections.^(1, 2) This organism is a significant cause of burn and nosocomial infections.⁽³⁾ The ability of *Pseudomonas aeruginosa* to destroy tissue may be related to the production of various extracellular enzymes.⁽²⁾

Pseudomonas aeruginosa produces a number of water-soluble pigments. Most of the strains of *Pseudomonas aeruginosa* produce pyocyanin, a blue, water- and chloroform-soluble, non-fluorescent pigment that diffuses into the surrounding medium. *Pseudomonas aeruginosa* is the only *Pseudomonas* species known to produce this pigment. (However, certain strains are apyocyanogenic).⁽⁴⁾ Some strains of *Pseudomonas aeruginosa* produce other pigments, such as the brown-black pyomelanin, the red pyorubin or the yellow pyoverdin. Pyoverdin is a water soluble fluorescent pigment often produced by *Pseudomonas aeruginosa* and other pseudomonads isolated from humans. When pyoverdin combines with the blue water-soluble pigment pyocyanin, the bright green color characteristic of *Pseudomonas aeruginosa* is created. The presence of these pigments can, however, mask the production of pyocyanin.

Pseudomonas Agar is based on the formulation described by King et al⁽⁵⁾ and as recommended in Indian Pharmacopoeia⁽⁶⁾ for detecting pyocyanin. This medium enhances the elaboration of pyocyanin but inhibits the formation of fluorescein pigment. Both pigments diffuse from *Pseudomonas* colonies into the medium when they grow. Pyocyanin elaborated on Pseudomonas Agar For Detection of Pyocyanin is a blue color.

Principles of the Procedure

Pseudomonas Agar For Detection of Pyocyanin contains pancreatic digest of gelatin which provides amino acids and other essential nitrogenous substances. Anhydrous potassium sulphate and magnesium chloride enhances the pyocyanin production and suppresses the fluorescein production. Low content of phosphorous in the medium also aids in inhibiting the production of fluorescein. Glycerol is supplemented as a source of carbon.

Formula / Liter

Ingredients	Gms / Liter
Pancreatic digest of gelatin	20.00
Anhydrous potassium sulphate	10.00
Anhydrous magnesium chloride	1.40
Agar	15.00
Final pH: 7.2 ± 0.2 at 25°C	
Formula may be adjusted and/or supplemented as required to meet performance specifications	

Precautions

1. For Laboratory Use only.
2. IRRITANT. Irritating to eyes, respiratory system, and skin.

Directions

1. Suspend 46.40 grams of the medium in one liter of distilled water containing 10 ml glycerol.
2. Heat to boiling, to dissolve the medium completely.
3. Autoclave at 121°C, 15 psi pressure, for 15 minutes / validated cycle.



PRODUCT SPECIFICATION SHEET

Quality Control Specifications

Dehydrated Appearance	Cream to yellow homogeneous free flowing powder
Prepared Medium	Yellow coloured clear to slightly opalescent gel forms in Petri plates
Reaction of % solution	Not applicable
Gel Strength	Firm, comparable with 1.5% Agar gel

Expected Cultural Response: Growth Promotion is carried out in accordance with the harmonized method of IP. Cultural response was observed after an incubation at 33-37°C for not less than 3 days. Recovery rate is considered as 100% for bacteria growth on Soyabean Casein Digest Agar.

Sr. No.	Organisms	Results to be achieved					Oxidase
		Inoculum (CFU)	Observed Lot value (CFU)	Recovery	Characteristic colonial morphology	Fluorescence in UV light	
	Test for <i>Pseudomonas aeruginosa</i>						
1.	<i>Pseudomonas aeruginosa</i> ATCC 9027	50-100	35 -100	≥70%	generally greenish	positive	positive
	Additional Microbiological Testing						
2.	<i>Pseudomonas aeruginosa</i> ATCC 27853	50-100	35 -100	≥70%	generally greenish	positive	positive

The organisms listed are the minimum that should be used for quality control testing.

Test Procedure

1. Obtain inoculum from a pure 18-24 hours culture of organism to be tested.
2. Inoculate plates or agar slants by streaking the surface.
3. Incubate at 35-37°C for 18-24 hours.
4. If the isolate fails to grow or grows slowly, reincubate at 25-30°C for 1-2 days and observe for growth and pigment production.

Results

1. A positive result is indicated by a blue pigment diffusing into the agar.
2. Pseudomonas Agar (For Pyocyanin), a blue to blue-green pigment seen in the colonies and surrounding medium. The pyocyanin pigment diffuses from the colonies of *Pseudomonas* into the agar and shows blue colouration.
3. Confirm the presence of pyocyanin by adding several drops of chloroform and observe for a blue color in the chloroform. (Pyocyanin is more soluble in chloroform than in water).

Storage

Store the sealed bottle containing the dehydrated medium at 10 - 30°C. Once opened and recapped, place container in a low humidity environment at the same storage temperature. Protect from moisture and light.

Expiration

Refer to the expiration date stamped on the container. The dehydrated medium should be discarded if not free flowing, or if the appearance has changed from the original color. Expiry applies to medium in its intact container when stored as directed.

Limitations of the Procedure

1. Occasionally, a *Pseudomonas* culture is encountered that will produce small amounts of pigment in the medium, indicated as a blue-green color on Pseudomonas Agar (For Pyocyanin). If a blue-green color occurs confirmation of the presence of pyocyanin can be made by extraction with chloroform (CHCl₃).
2. The formation of nonpigmented colonies does not completely rule-out a *Pseudomonas aeruginosa* isolate.
3. A pyocyanin-producing *Pseudomonas* strain will usually also produce fluorescein. It must, therefore, be differentiated from other simple fluorescent pseudomonads by other means. Temperature can be a determining factor as most other fluorescent strains will not grow at 35°C. Rather, they grow at 25-30°C.

Packaging

Product Name : Pseudomonas Agar Medium For Detection of Pyocyanin

Product Code : DM219I

Available Pack sizes : 100gm/ 500gm



PRODUCT SPECIFICATION SHEET

References

1. Kiska and Gilligan. 1999. *In* Murray, Baron, Pfaller, Tenover and Tenover (ed.), Manual of clinical microbiology, 7th ed. American Society for Microbiology, Washington, D.C.
2. Baron, E. J., L. R. Peterson, and S. M. Finegold. 1994. Nonfermentative gram-negative bacilli and coccobacilli, p. 386-405. Bailey & Scott's diagnostic microbiology, 9th ed. Mosby-Year Book, Inc. St. Louis, MO.
3. Gilligan, P. H. 1995. *Pseudomonas* and *Burkholderia*, p. 509-519. *In* P. R. Murray, E. J. Baron, M. A. Pfaller, F. C. Tenover, and R. H. Tenover (eds.), Manual of clinical microbiology, 6th ed. American Society of Microbiology, Washington, D.C.
4. Forbes, Sahn and Weissfeld. 2007. Bailey & Scott's diagnostic microbiology, 12th ed. Mosby, Inc., St. Louis, Mo.
5. King, Ward and Raney, 1954, J.Lab. and Clin. Med., 44:301
6. Indian Pharmacopoeia, 2007, Govt. of India, Ministry of Health and Family Welfare, New Delhi, India.

Further Information

For further information please contact your local MICROMASTER Representative.



MICROMASTER LABORATORIES PRIVATE LIMITED

DM219IPSS, Rev.00, Ver.00/ 01.02.2016

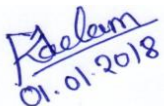
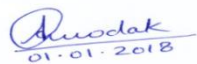

Unit 38/39, Kalpataru Industrial Estate,

Off G.B. Road, Near 'R-Mall', Thane (W) - 400607. M.S. INDIA.

Ph: +91-22-25895505, 4760, 4681. Cell: 9320126789.

Email: micromaster@micromasterlab.com

sales@micromasterlab.com

	Checked By	Approved By
 01.01.2018	 01.01.2018	 01.01.2018
Microbiologist	Head Quality Control	Head Quality Assurance

Disclaimer :

All Products conform exclusively to the information contained in this and other related Micromaster Publications. Users must ensure that the product(s) is appropriate for their application, prior to use. The information published in this publication is based on research and development work carried out in our laboratory and is to the best of our knowledge true and accurate. Micromaster Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are intended for laboratory, diagnostic, research or further manufacturing use only and not for human or animal or therapeutic use, unless otherwise specified. Statements included herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.