



PRODUCT SPECIFICATION SHEET

Simmons Citrate Agar (DM244)

Intended Use

Simmons Citrate Agar (DM244) is used for differentiating members of Enterobacteriaceae on the basis of citrate utilization.

Product Summary and Explanation

Simmons Citrate Agar is recommended by Ewing and Edwards⁽¹⁾ for the differentiation of the family Enterobacteriaceae based on whether or not citrate is utilized as the sole source of carbon. Koser⁽²⁾ first developed a liquid medium for differentiating coliforms. Fecal coliforms were unable to use citrate as the sole source of carbon and inorganic ammonium salt as a sole source of nitrogen. Non-fecal coliform, such as *Enterobacter aerogenes* or *Salmonella enteritidis* could use citrate in such a medium with resultant alkalinity. Simmons Citrate Agar is a modification of Koser's medium to which bromothymol blue and 1.5% agar have been added. Organisms able to metabolize the citrate grow luxuriantly. The medium is alkalized and changes from its initial green to deep blue in 24-48 hours. E. coli either do not grow at all on this medium, or grow so sparsely that no change in reaction is apparent. Simmons citrate agar is recommended for differentiation of enteric gram-negative bacilli from clinical specimens,^(3,4) water samples,⁽⁵⁾ and food samples.^(6,7,8,9)

Principles of the Procedure

Ammonium Dihydrogen Phosphate is the sole source of nitrogen in Simmons Citrate Agar. Magnesium is a cofactor for a variety of metabolic reactions. Phosphate acts as a buffer. Sodium citrate is the sole source of carbon in this medium. Sodium chloride maintains the osmotic balance of the medium. Agar is the solidifying agent. Bromothymol blue is the pH indicator. Organisms that can utilize ammonium dihydrogen phosphate and sodium citrate as their sole sources of nitrogen and carbon will grow on this medium and produce a color change from green (neutral) to blue (alkaline).

Formula / Liter

Ingredients	Gms / Liter
Magnesium sulphate	0.20
Ammonium dihydrogen phosphate	1.00
Dipotassium phosphate	1.00
Sodium citrate	2.00
Sodium chloride	5.00
Bromothymol blue	0.08
Agar	15.00
Final pH (at 25°C) 6.8 ± 0.2	
Formula may be adjusted and/or supplemented as required to meet performance specifications	

Precautions

1. For Laboratory Use only.
2. IRRITANT. Irritating to eyes, respiratory system, and skin.

Directions

1. Suspend 24.28grams of the medium in one liter of distilled water.
2. Heat if necessary, to dissolve the medium completely.
3. Autoclave at 121°C , 15 psi pressure, for 15 minutes / validated cycle.
4. Mix well and pour into sterile Petri plates.





PRODUCT SPECIFICATION SHEET


Quality Control Specifications

Dehydrated Appearance	Cream to yellow colored, homogeneous, free flowing powder
Prepared Medium	Forest green coloured slightly opalescent gel forms in tubes as slants
Reaction of 2.43% Solution	pH 6.8 ± 0.2 at 25°C
Gel Strength	Firm, comparable with 1.5% Agar gel

Expected Cultural Response: Cultural characteristics observed after an incubation at 35 - 37°C for 18 - 48 hours.

Sr. No.	Organisms	Results to be achieved		
		Inoculum (CFU)	Growth	Citrate utilisation
1.	<i>Enterobacter aerogenes</i> ATCC 13048	50-100	Good-luxuriant	Positive reaction, blue colour
2.	<i>Escherichia coli</i> ATCC 25922	>=10 ³	Inhibited	--
3.	<i>Salmonella choleraesuis</i> ATCC 12011	50-100	Good-luxuriant	Positive reaction, blue colour
4.	<i>Salmonella enteritidis</i> ATCC 13076	50-100	Good-luxuriant	Positive reaction, blue colour
5.	<i>Salmonella typhi</i> ATCC 6539	50-100	Fair-good	Negative reaction, green colour
6.	<i>Salmonella typhimurium</i> ATCC 14028	50-100	Good-luxuriant	Positive reaction, blue colour
7.	<i>Shigella dysenteriae</i> ATCC 13313	>=10 ³	Inhibited	--

The organisms listed are the minimum that should be used for quality control testing.

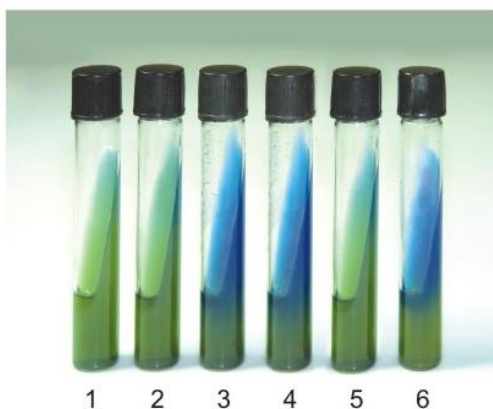

Simmons Citrate Agar (DM244)

1. Control
2. *Escherichia coli* ATCC 25922
3. *Enterobacter aerogenes* ATCC 13048
4. *Salmonella typhimurium* ATCC 14028
5. *Salmonella typhi* ATCC 6539
6. *Salmonella enteritidis* ATCC 13076

Test Procedure

1. Obtain a pure culture of the organism to be tested.
2. With an inoculating needle or loop, pick the center of a well-isolated colonies obtained from solid culture media.
3. Streak only the surface of the slant with a light inoculum.
4. Incubate at 35± 2°C for 18-48 hours.





PRODUCT SPECIFICATION SHEET

Results

A Positive reaction is indicated by growth on the slant with an intense blue color (alkaline reaction)

Storage

Store the sealed bottle containing the dehydrated medium at 2 - 30°C. Once opened and recapped, place container in a low humidity environment at the same storage temperature. Protect from moisture and light.

Expiration

Refer to the expiration date stamped on the container. The dehydrated medium should be discarded if not free flowing, or if the appearance has changed from the original color. Expiry applies to medium in its intact container when stored as directed.

Limitations of the Procedure

1. Due to nutritional variation, some strains may be encountered that grow poorly or fail to grow on this medium.
2. Some citrate positive organisms require a 48 hour incubation or longer for a pH change to occur.⁽¹⁰⁾

Packaging

Product Name: Simmons Citrate Agar

Product Code: DM244

Available Pack sizes: 100gm / 500gm

References

1. Edwing W. H. and Edwards P.R. (1960) Bull. Bact. Nomen. and Taxon. 10. 1-12.
2. Koser, S.A. 1923. Utilization of the salts of organic acids by the colon-aerogenes group. J. Bacteriol 8:493.
3. Pezzlo, M. (ed). 1992. Aerobic bacteriology, p. 1.0.0-1.20.47. In Isenberg, H.D. (ed), Clinical microbiology procedures handbook, vol. 1. American Society for Microbiology, Washington, D.C.
4. Baron, E. J., L. R. Peterson, S.M. Finegold. 1994. Bailey & Scott's diagnostic microbiology, 9th ed. Mosby-Year Book, Inc., St. Louis, MO.
5. Eaton, A.D., L.S. Clesceri, and A.E. Greenberg (ed.). 1995. Standard methods for the examination of water and wastewater, 19th ed. American Public Health Association, Washington, D.C.
6. Vanderzant, C., and D.F. Splittstoesser (ed.). 1992. Compendium of methods for the microbiological examination of foods, 3rd ed. American Public Health Association, Washington, D.C.
7. FDA Bacteriological Analytical Manual, 8th ed. AOAC International, Gaithersburg, MD.
8. Association of Official Analytical Chemists. 1995. Official methods of analysis of AOAC International, 16th ed. AOAC International, Arlington, VA.
9. Federal Register. 1996. Pathogen reduction; hazard analysis and critical point (HACCP) systems; final rule. Fed. Regis. 61:38917-38925.
10. MacFaddin, J.D. 1985. Media for isolation -cultivation-identification-maintenance of medical bacteria, vol. 1. Williams & Wilkins, Baltimore, MD.

Further Information

For further information please contact your local MICROMASTER Representative.





PRODUCT SPECIFICATION SHEET






MICROMASTER LABORATORIES PRIVATE LIMITED

DM244PSS,QAD/FR/024,Rev.00/01.01.2018

Unit 38/39, Kalpataru Industrial Estate,
Off G.B. Road, Near 'R-Mall', Thane (W) - 400607. M.S. INDIA.
Ph: +91-22-25895505, 4760, 4681. Cell: 9320126789.

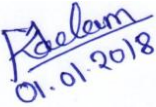Email: micromaster@micromasterlab.com
sales@micromasterlab.com

Prepared By	Checked By	Approved By
 01.01.2018	 01.01.2018	 01.01.2018
Microbiologist	Head Quality Control	Head Quality Assurance

Disclaimer :

All Products conform exclusively to the information contained in this and other related Micromaster Publications. Users must ensure that the product(s) is appropriate for their application, prior to use. The information published in this publication is based on research and development work carried out in our laboratory and is to the best of our knowledge true and accurate. Micromaster Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are intended for laboratory, diagnostic, research or further manufacturing use only and not for human or animal or therapeutic use, unless otherwise specified. Statements included herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.

