



## PRODUCT SPECIFICATION SHEET

### Casein Soyabean Digest Agar (Soyabean Casein Digest Agar) (DM247B)

#### Intended Use

Casein Soyabean Digest Agar (Soyabean Casein Digest Agar) (DM247B) is recommended for cultivation of a wide variety of microorganisms from pharmaceutical products, in compliance with BP.

#### Product Summary and Explanation

Soyabean Casein Digest Agar is recommended as sterility testing medium by various pharmacopoeias. This medium is prepared as per British Pharmacopeia.<sup>(2)</sup> This medium is also employed in validation of sterility checking procedure in accordance with the microbial limit testing harmonized methodology of USP/BP/EP/JP/IP.<sup>(1,2,3,4,5)</sup> This medium is used in microbial limit test and antimicrobial preservative- effective test. Gunn et al<sup>(6)</sup> applied this medium for the growth of fastidious organisms and for study of haemolytic reaction after addition of 5% v/v blood.

#### Principles of the Procedure

Soyabean Casein Digest Agar contains pancreatic digest of casein and papaic digest of soyabean meal which makes these media nutritious by providing amino acids and long chain peptides for the growth of microorganisms. Natural sugars of soy enhance growth of microorganism. Sodium chloride maintains the osmotic balance. Agar is the solidifying agent.

#### Formula / Liter

Ingredients	Gms / Liter
Pancreatic digest of casein	15.00
Papaic digest of soyabean meal	5.00
Sodium chloride	5.00
Agar	15.00
Final pH: 7.3 ± 0.2 at 25°C	
Formula may be adjusted and/or supplemented as required to meet performance specifications	

#### Precautions

1. For Laboratory Use only.
2. IRRITANT. Irritating to eyes, respiratory system, and skin.

#### Directions

1. Suspend 40 grams of the medium in one liter of distilled water.
2. Heat to boiling, to dissolve the medium completely.
3. Autoclave at 121°C, 15 psi pressure, for 15 minutes / validated cycle.
4. Mix well and pour into sterile petri plates.

#### Quality Control Specifications

Dehydrated Appearance	Cream to yellow homogeneous free flowing powder
Prepared Medium	Light yellow coloured clear to slightly opalescent gel forms in Petri plates
Reaction of % Solution	Not Applicable
Gel Strength	Firm, comparable with 1.5% Agar gel





## PRODUCT SPECIFICATION SHEET

### Growth Promotion Test

Growth Promotion was carried out in accordance with the harmonized method of BP, and growth was observed after an incubation at 30-35°C for 18-24 hours. Recovery rate is considered 100% for bacteria growth on Blood Agar and fungus growth on Sabouraud Dextrose Agar.

### Growth promoting properties

Growth of microorganism comparable to that previously obtained with previously tested and approved lot of medium occurs at the specified temperature for not more than the shortest period of time specified inoculating  $\leq 100$  cfu (at 30-35°C for 18 hours).

### Expected Cultural Response:

Sr. No.	Organisms	Results to be achieved CFU		
		Inoculum (CFU)	Observed Lot value (CFU)	Recovery
1.	<i>Bacillus subtilis</i> ATCC 6633	50 - 100	35 -100	$\geq 70$ %
2.	<i>Staphylococcus aureus</i> ATCC 25923	50 - 100	35 -100	$\geq 70$ %
3.	<i>Staphylococcus aureus</i> ATCC 6538	50 - 100	35 -100	$\geq 70$ %
4.	<i>Escherichia coli</i> ATCC 25922	50 - 100	35 -100	$\geq 70$ %
5.	<i>Escherichia coli</i> ATCC 8739	50 - 100	35 -100	$\geq 70$ %
6.	<i>Escherichia coli</i> NCTC 9002	50 - 100	35 -100	$\geq 70$ %
7.	<i>Pseudomonas aeruginosa</i> ATCC 27853	50 - 100	35 -100	$\geq 70$ %
8.	<i>Pseudomonas aeruginosa</i> ATCC 9027	50 - 100	35 -100	$\geq 70$ %
9.	<i>Salmonella</i> Abony NCTC 6017	50 - 100	35 -100	$\geq 70$ %
10.	<i>Micrococcus luteus</i> ATCC 9341	50 - 100	35 -100	$\geq 70$ %
11.	<i>Streptococcus pneumonia</i> ATCC 6305	50 - 100	35 -100	$\geq 70$ %
12.	<i>Salmonella</i> Typhimurium ATCC 14028	50 - 100	35 -100	$\geq 70$ %
13.	<i>Candida albicans</i> ATCC 10231	50 - 100	35 -100	$\geq 70$ %
14.	<i>Candida albicans</i> ATCC 2091	50 - 100	35 -100	$\geq 70$ %
15.	<i>Aspergillus brasiliensis</i> ATCC 16404	50 - 100	25 -70	50-70%

The organisms listed are the minimum that should be used for quality control testing.

### Test Procedure

Refer to appropriate references for standard test procedures.

### Results

1. After incubation, it is desirable to have isolated colonies of organisms from the original sample. Subculture colonies of interest, so that positive identification can be made by means of biochemical and/or serological testing.





## PRODUCT SPECIFICATION SHEET

---

2. The total aerobic count is considered to be equal to the number of colony forming units found on this medium, if colonies of fungi are detected on this medium they are counted along with total aerobic count.

### Storage

Store the sealed bottle containing the dehydrated medium at 10 - 30°C. Once opened and recapped, place container in a low humidity environment at the same storage temperature. Protect from moisture and light.

### Expiration

Refer to the expiration date stamped on the container. The dehydrated medium should be discarded if not free flowing, or if the appearance has changed from the original color. Expiry applies to medium in its intact container when stored as directed.

### Limitations of the Procedure

1. For identification, organisms must be in pure culture. Morphological, biochemical and/or serological tests should be performed for final identification.
2. Consult appropriate texts for detailed information and recommended procedures.

### Packaging

**Product Name : Casein Soyabean Digest Agar (Soyabean Casein Digest Agar)**

**Product Code : DM247B**

**Available Pack sizes : 100gm/ 500gm**

### References

1. The United States Pharmacopoeia, 2011, The United States Pharmacopoeial Convention. Rockville, MD.
2. British Pharmacopoeia, 2011, The Stationery office British Pharmacopoeia
3. European Pharmacopoeia, 2011, European Dept. for the quality of Medicines.
4. Japanese Pharmacopoeia, 2008.
5. The Indian Pharmacopoeia 2010, Govt of India, Ministry of Health and Family Welfare, New Delhi.
6. Gunn. B. A. et al, 1977, J. Clin. Microbiol., 5(6) : 650

### Further Information

For further information please contact your local MICROMASTER Representative.



**MICROMASTER LABORATORIES PRIVATE LIMITED**

DM247BPSS,QAD/FR/024,Rev.00/01.01.2018

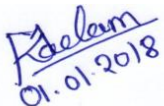

Unit 38/39, Kalpataru Industrial Estate,  
Off G.B. Road, Near 'R-Mall' , Thane (W) - 400607. M.S. INDIA.  
Ph: +91-22-25895505, 4760, 4681. Cell: 9320126789.  
Email: [micromaster@micromasterlab.com](mailto:micromaster@micromasterlab.com)  
[sales@micromasterlab.com](mailto:sales@micromasterlab.com)





## PRODUCT SPECIFICATION SHEET

---

Prepared By	Checked By	Approved By
		
Microbiologist	Head Quality Control	Head Quality Assurance

**Disclaimer :**

All Products conform exclusively to the information contained in this and other related Micromaster Publications. Users must ensure that the product(s) is appropriate for their application, prior to use. The information published in this publication is based on research and development work carried out in our laboratory and is to the best of our knowledge true and accurate. Micromaster Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are intended for laboratory, diagnostic, research or further manufacturing use only and not for human or animal or therapeutic use, unless otherwise specified. Statements included herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.

