



## PRODUCT SPECIFICATION SHEET

### Bile Salt Agar (DM299BS)

#### Intended Use

Bile Salt Agar (DM299BS) is recommended for isolation and enumeration of food poisoning-bile tolerant enteric bacilli, in compliance with Indian Standard published by BIS. IS:5887 (Part V)-1976.

#### Product Summary and Explanation

*V. cholerae* thrives in water ecology, particularly surface water. Cholera infections are most commonly acquired from drinking water in which *V. cholerae* is found naturally or into which it has been introduced from the feces of an infected person. *Vibrio* species, like many other gram-negative bacteria, grow in the presence of relatively high levels of bile salts.<sup>(1)</sup> Bile Salt Agar is formulated as recommended by BIS<sup>(2)</sup> for isolation, identification and enumeration of *Vibrio cholerae*. On Bile Salt Agar colonies of *Vibrio* species have a distinctive appearance which may be seen by growing a known strain of *Vibrio cholerae* and comparing it with *Escherichia coli*.<sup>(2)</sup> Slide agglutination test is carried out using polyvalent cholera typing serum for suspicious growths.

#### Principles of the Procedure

Bile Salt Agar contains peptic digest of animal tissue, meat extract which provides carbonaceous, nitrogenous compounds and other essential nutrients growth for enteric bacilli. Sodium chloride maintains the osmotic balance of the medium. Sodium taurocholate inhibits most of the gram-negative organisms.

#### Formula / Liter

Ingredients	Gms / Liter
Peptic digest of animal tissue	10.00
Meat extract	5.00
Sodium chloride	5.00
Sodium taurocholate	5.00
Agar	15.00
Final pH: 8.5 ± 0.2 at 25°C	
Formula may be adjusted and/or supplemented as required to meet performance specifications	

#### Precautions

1. For Laboratory Use only.
2. IRRITANT. Irritating to eyes, respiratory system, and skin.

#### Directions

1. Suspend 40 grams of the medium in one liter of distilled water.
2. Heat if necessary, to dissolve the medium completely.
3. Autoclave at 121°C, 15 psi pressure, for 15 minutes / validated cycle.
4. Mix well and pour into sterile petri plates.





## PRODUCT SPECIFICATION SHEET

### Quality Control Specifications

Dehydrated Appearance	Cream to yellow coloured homogeneous free flowing powder
Prepared Medium	Light yellow coloured clear to slightly opalescent gel forms in petri plates
Reaction of 4.0% Solution	pH : 8.5 ± 0.2 at 25°C
Gel Strength	Firm, comparable with 1.5% Agar gel

**Expected Cultural Response:** Cultural characteristics observed after an incubation at 35-37°C for 18-24 hours.

Sr. No.	Organisms	Results to be achieved	
		Inoculum (CFU)	Growth
1.	<i>Enterobacter aerogenes</i> ATCC 13048	50-100	good-luxuriant
2.	<i>Escherichia coli</i> ATCC 25922	50-100	good-luxuriant
3.	<i>Salmonella Typhi</i> ATCC 6539	50-100	good-luxuriant
4.	<i>Staphylococcus aureus</i> ATCC 25923	50-100	good-luxuriant
5.	<i>Vibrio cholerae</i> ATCC 15748	50-100	good-luxuriant

The organisms listed are the minimum that should be used for quality control testing.

### Test Procedure

Refer appropriate references for specific test procedures for isolation and enumeration of enteric bacilli.

### Results

Refer appropriate references and test procedures for interpretation of results.

### Storage

Store the sealed bottle containing the dehydrated medium at 10 - 30°C. Once opened and recapped, place container in a low humidity environment at the same storage temperature. Protect from moisture and light.

### Expiration

Refer to the expiration date stamped on the container. The dehydrated medium should be discarded if not free flowing, or if the appearance has changed from the original color. Expiry applies to medium in its intact container when stored as directed.

### Limitations of the Procedure

1. For identification, organisms must be in pure culture. Morphological, biochemical and/or serological tests should be performed for final identification.
2. Consult appropriate texts for detailed information and recommended procedures.

### Packaging

Product Name : Bile Salt Agar

Product Code : DM299BS

Available Pack sizes : 100gm / 500gm





## PRODUCT SPECIFICATION SHEET

### References

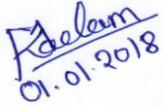

1. Vanderzant C. and Splittstoesser D. (Eds.), 1992, Compendium of Methods for the Microbiological Examination of Foods, 3rd ed., APHA, Washington, DC.
2. Bureau of Indian Standards IS : 5887 (Part V) Reaffirmed 1986.

### Further Information

For further information please contact your local MICROMASTER Representative.



**MICROMASTER LABORATORIES PRIVATE LIMITED**      DM299BSPSS,QAD/FR/024,Rev.00/01.01.2018  
Unit 38/39, Kalpataru Industrial Estate,  
Off G.B. Road, Near 'R-Mall', Thane (W) - 400607. M.S. INDIA.  
Ph: +91-22-25895505, 4760, 4681. Cell: 9320126789.  
Email: [micromaster@micromasterlab.com](mailto:micromaster@micromasterlab.com)  
[sales@micromasterlab.com](mailto:sales@micromasterlab.com)

Prepared By	Checked By	Approved By
 01.01.2018	 01.01.2018	 01.01.2018
Microbiologist	Head Quality Control	Head Quality Assurance





## PRODUCT SPECIFICATION SHEET

---

### Disclaimer :

All Products conform exclusively to the information contained in this and other related Micromaster Publications. Users must ensure that the product(s) is appropriate for their application, prior to use. The information published in this publication is based on research and development work carried out in our laboratory and is to the best of our knowledge true and accurate. Micromaster Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are intended for laboratory, diagnostic, research or further manufacturing use only and not for human or animal or therapeutic use, unless otherwise specified. Statements included herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.

