



## PRODUCT SPECIFICATION SHEET

### Rose Bengal Agar Base (DM327)

#### Intended Use

Rose Bengal Agar Base (DM326) is recommended for selective isolation and enumeration of yeast and mould from environmental material and food stuff.

#### Product Summary and Explanation

Various methods have been described for the selective isolation of fungi from environmental materials and foodstuffs containing mixed populations of fungi and bacteria. The use of media with an acid pH that selectively inhibits the growth of bacteria and thereby promotes the growth of fungi has been widely employed.<sup>(1-3)</sup> A number of investigators have reported, that, neutral pH media with antibiotics is advantageous for fungal growth compared to acidified media as the later may actually inhibit fungal growth,<sup>(4,5)</sup> fail to completely inhibit bacterial growth<sup>(5)</sup> and have little effect in restricting the size of mold colonies.<sup>(6)</sup> Smith and Dawson<sup>(7)</sup> used rose bengal in a neutral pH medium for the selective isolation of fungi from soil samples. Chloramphenicol, streptomycin, oxytetracycline and chlortetracycline have been used for the improved, selective isolation and enumeration of yeasts and molds from soil, sewage and foodstuffs.<sup>(4, 8-11)</sup> Rose Bengal Agar Base supplemented with Rose Bengal Chloramphenicol is a modification of the Rose Bengal Chlortetracycline Agar formula of Jarvis.<sup>(11)</sup> Chloramphenicol is employed instead of chlortetracycline, in this medium as a selective supplement. Of the antibiotics most frequently employed in media of neutral pH, chloramphenicol is recommended because of its heat stability and broad antibacterial spectrum. A modified formulation of Rose Bengal Agar is recommended in *Standard Methods for the Examination of Water and Wastewater* for the enumeration of yeasts and molds.<sup>(12)</sup>

#### Principles of the Procedure

Rose Bengal Agar Base contains papaic digest of soyabean meal which provides the carbon and nitrogen sources required for good growth of a wide variety of organisms. Dextrose is an energy source. Monopotassium phosphate provides buffering capability. Magnesium sulphate provides necessary trace elements. Rose bengal is a selective agent that inhibits bacterial growth and restricts the size and height of colonies of the more rapidly growing molds. The restriction in growth of molds aids in the isolation of slow-growing fungi by preventing overgrowth by more rapidly growing species. Rose bengal is taken up by yeast and mold colonies, thereby facilitating their recognition and enumeration. Chloramphenicol Selective Supplement (MS057) inhibit bacteria.

#### Formula / Liter

Ingredients	Gms / Liter
Papaic digest of soyabean meal	5.00
Dextrose	10.00
Monopotassium phosphate	1.00
Magnesium sulphate	0.50
Rose Bengal	0.05
Agar	15.00
Final pH: 7.2 ± 0.2 at 25°C	
Formula may be adjusted and/or supplemented as required to meet performance specifications	

#### Precautions

1. For Laboratory Use only.
2. IRRITANT. Irritating to eyes, respiratory system, and skin.

#### Directions

1. Suspend 31.55 grams of the medium in one liter of distilled water.



## PRODUCT SPECIFICATION SHEET

- Heat to boiling, to dissolve the medium completely.
- Autoclave at 121°C, 15 psi pressure, for 15 minutes / validated cycle.
- Cool to 40-45°C and add 2 ml of rehydrated Chloramphenicol Selective Supplement (MS057) for each 500 ml of medium.
- Mix thoroughly and pour into sterile Petri plates.

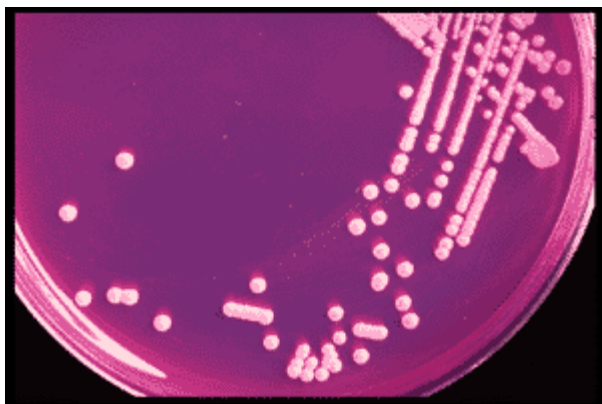
### Quality Control Specifications

Dehydrated Appearance	Light yellow to pink homogeneous free flowing powder
Prepared Medium	Deep pink coloured clear to very slightly opalescent gel forms in Petri plates
Reaction of 3.16% Solution	pH : 7.2 ± 0.2 at 25°C
Gel Strength	Firm, comparable with 1.5% Agar gel

**Expected Cultural Response:** Cultural characteristics observed after an incubation at 20-25°C for 5 days with added Chloramphenicol Selective Supplement (MS057).

Sr. No.	Organisms	Results to be achieved		
		Inoculum (CFU)	Growth	Recovery
1.	<i>Aspergillus brasiliensis</i> ATCC 16404	50 -100	good	--
2.	<i>Candida albicans</i> ATCC 10231	50 -100	good	≥50%
3.	<i>Escherichia coli</i> ATCC 25922	≥10 <sup>3</sup>	inhibited	0%
4.	<i>Micrococcus luteus</i> ATCC 10240	≥10 <sup>3</sup>	inhibited	0%
5.	<i>Saccharomyces cerevisiae</i> ATCC 9763	50 -100	good	≥50%

The organisms listed are the minimum that should be used for quality control testing.



**Rose Bengal Agar (DM327)**  
*Candida albicans* ATCC 10231

### Test Procedure

- Pour plates : pipette 1 ml aliquots of a suitable series of dilution to Petri plates. Add 10-12 ml of molten medium, cooled at 45°C, and mix thoroughly.  
Spread plates: spread 0.1ml of aliquots of a suitable series of dilution over the surface of the solidified medium in a Petri dish.
- Incubate for upto 5 days at 25°C.



## PRODUCT SPECIFICATION SHEET

### Results

1. Calculate the number of yeasts or moulds per 1 gram or 1 ml by multiplying the number of colonies by dilution factor.
2. Colonies of yeast appear pink due to uptake of rose bengal.

### Storage

Store the sealed bottle containing the dehydrated medium at 10 - 30°C. Once opened and recapped, place container in a low humidity environment at the same storage temperature. Protect from moisture and light.

### Expiration

Refer to the expiration date stamped on the container. The dehydrated medium should be discarded if not free flowing, or if the appearance has changed from the original color. Expiry applies to medium in its intact container when stored as directed.

### Limitations of the Procedure

1. Although this medium is selective primarily for fungi, microscopic examination is recommended for presumptive identification. Biochemical testing using pure cultures is required for complete identification.
2. Due to the selective properties of this medium and the type of specimen being cultured, some strains of fungi may be encountered that fail to grow or grow poorly on the complete medium; similarly, some strains of bacteria may be encountered that are not inhibited or only partially inhibited.
3. Care should be taken not to expose this medium to light, since photo degradation of rose bengal yields compounds that are toxic to fungi.

### Packaging

**Product Name : Rose Bengal Agar Base**

**Product Code : DM327**

**Available Pack sizes : 500gm**

### References

1. Waksman S. A., 1922, J. Bacteriol., 7:339.
2. Koburger J. A., 1976, In Speck M. L., (Ed.), Compendium of Methods for the Microbiological Examination of Foods, American Public Health Association, Washington, D.C.
3. Mossel D. A. A., Visser M. and Mengerink W. H. J., 1962, Lab Practice 11:109.
4. Martin J. P., 1950, Soil Sci. 69:215.
5. Koburger J. A., 1972, J. Milk Food Technol. 35:659.
6. Tyner L. E., 1944, Soil Sci. 57:271.
7. Smith M. R. and Dawson V. T., 1944, Soil Sci. 58:467.
8. Cooke W. B., 1954, Antibiot. and Chemother., 4:657.
9. Papavizas and Davey, 1959. Soil Sci. 88:112.
10. Overcast W. W. and Weakley D. J., 1969, J. Milk Technol., 32:442.
11. Jarvis B., 1973, J. Appl. Bacteriol., 36:723.
12. Eaton, Rice and Baird (ed.). 2005. Standard methods for the examination of water and wastewater, 21<sup>st</sup> ed., online. American Public



**MICROMASTER LABORATORIES PRIVATE LIMITED**

Unit 38/39, Kalpataru Industrial Estate,  
Near Runwal Estate, Behind 'R-Mall', Ghodbunder Raod,  
Thane (W) - 400607. M.S. INDIA.

Ph: +91-22-25895505, 4760, Cell: 9320126789.

Email: [micromaster@micromasterlab.com](mailto:micromaster@micromasterlab.com)

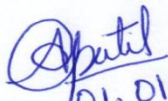
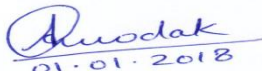

[sales@micromasterlab.com](mailto:sales@micromasterlab.com)

DM327PSSQAD/FR/024, Rev.00/01.01.2018





## PRODUCT SPECIFICATION SHEET

Prepared By	Checked By	Approved By
 01.01.2018	 01.01.2018	 01.01.2018
Microbiologist	Head Quality Control	Head Quality Assurance

### Disclaimer :

All Products conform exclusively to the information contained in this and other related Micromaster Publications. Users must ensure that the product(s) is appropriate for their application, prior to use. The information published in this publication is based on research and development work carried out in our laboratory and is to the best of our knowledge true and accurate. Micromaster Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are intended for laboratory, diagnostic, research or further manufacturing use only and not for human or animal or therapeutic use, unless otherwise specified. Statements included herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.

