



PRODUCT SPECIFICATION SHEET

Selenite Broth Base w/o Selenite (DM402)

Intended Use

Selenite Broth Base w/o Selenite (DM402) is recommended for enrichment of *Salmonellae* species from food, dairy products, pathological materials, with addition of selenite.

Product Summary and Explanation

Selective inhibitory effects of selenite were first demonstrated by Klett.⁽¹⁾ Guth⁽²⁾ used it to isolate *Salmonella typhi*. Leifson⁽³⁾ found that selenite inhibited fecal streptococci and coliforms during the first 8-12 hours of incubation, thereby permitting salmonellae to replicate without overwhelming interference from other members of the intestinal flora. Selenite Broth is useful for detecting *Salmonella* in the non-acute stages of illness when organisms occur in the faeces in low numbers and for epidemiological studies to enhance the detection of low numbers of organisms from asymptomatic or convalescent patients.⁽⁴⁾ *Salmonella* are also injured during various food processing procedures, including exposure to low temperatures, submarginal heat, drying, radiation, preservatives or sanitizers.⁽⁵⁾ Since, *Salmonella* may be present in low numbers in food sample in an injured condition; recovery of *Salmonella* involves pre-enrichment, selective enrichment and selective plating..

Principles of the Procedure

Selenite Broth Base w/o Selenite contains casein enzymic hydrolysate provides nitrogenous substances. Lactose is the fermentable carbohydrate and maintains the pH in medium as selenite is reduced by bacterial growth and alkali is produced. An increase in pH lowers the toxicity of the selenite and results in overgrowth of other bacteria. The acid produced by bacteria due to lactose fermentation serves to maintain a neutral pH. Enriched broth is subcultured on solid medium. Do not incubate the broth longer than 24 hours as inhibitory effect of selenite reduces after 6 - 12 hours of incubation.⁽⁶⁾

Formula / Liter

Ingredients	Gms / Liter
Casein enzymic hydrolysate	5.00
Lactose	4.00
Sodium phosphate	10.00
Final pH: 7.0 ± 0.2 at 25°C	
Formula may be adjusted and/or supplemented as required to meet performance specifications	

Precautions

1. For Laboratory Use only.
2. IRRITANT. Irritating to eyes, respiratory system, and skin.
3. Sodium hydrogen selenite (Sodium bi-selenite) is very toxic, corrosive agent and causes teratogenicity. Handle with great care. Upon contact with skin, wash immediately with a lot of water.

Directions

1. Suspend 4 grams of Sodium Biselenite(MS050) in one liter of distilled water.
2. Add 19.01 grams of medium. Mix well.
3. Warm to dissolve the medium completely.
4. Distribute in sterile test tubes.





PRODUCT SPECIFICATION SHEET

5. Sterilize in a boiling water bath or free flowing steam for 10 minutes.
6. DO NOT AUTOCLAVE.
7. Excessive heating is detrimental. Discard the prepared medium if large amount of selenite is reduced (indicated by red precipitate at the bottom of tube / bottle).

Quality Control Specifications

Dehydrated Appearance	Part A :Cream to Light yellow homogeneous free flowing powder
Prepared Medium	Light yellow coloured, clear solution without any precipitate
Reaction of [(1.9% w/v) Medium and (0.4% w/v) Sodium Biselenite] Solution	pH : 7.0 ± 0.2 at 25°C
Gel Strength	Not Applicable

Expected Cultural Response: Cultural characteristics observed after an incubation at 35-37°C for 18-24 hours when sub cultured on MacConkey Agar.

Sr. No.	Organisms	Results to be achieved		
		Inoculum (CFU)	Growth	Colour of the medium
1.	<i>Escherichia coli</i> ATCC 25922	50-100	little-none (no increase in numbers)	pink with bile precipitate
2.	<i>Salmonella choleraesuis</i> ATCC 12011	50-100	good-luxuriant	colourless
3.	<i>Salmonella typhimurium</i> ATCC 14028	50-100	good-luxuriant	colourless
4.	<i>Salmonella typhi</i> ATCC 6539	50-100	good-luxuriant	colourless

The organisms listed are the minimum that should be used for quality control testing.

Test Procedure

1. Prepare food sample following the recommended procedure.
2. Inoculate into recommended pre-enrichment broth.
3. Transfer 1 mL of mixture to 10 mL Selenite Cystine Broth and to 10 mL Tetrathionate Broth.
4. Incubate at 35°C for 24 ± 2 hours.
5. Mix and streak 3 mm loopful (10 µL) of sample from both broths onto Bismuth Sulfite Agar, Xylose Lysine Desoxycholate Agar, Hektoen Enteric Agar or MacConkey Agar.
6. Examine plates for the presence of colonies that are typical for *Salmonella* spp.
7. Refer appropriate references for standard test procedures.

Results

Refer appropriate references and procedures for interpretation of results.

Storage





PRODUCT SPECIFICATION SHEET

Store the sealed bottle containing the dehydrated medium at 10 - 30°C. Once opened and recapped, place container in a low humidity environment at the same storage temperature. Protect from moisture and light.

Expiration

Refer to the expiration date stamped on the container. The dehydrated medium should be discarded if not free flowing, or if the appearance has changed from the original color. Expiry applies to medium in its intact container when stored as directed.

Limitations of the Procedure

1. For identification, organisms must be in pure culture. Morphological, biochemical and/or serological tests should be performed for final identification.
2. Consult appropriate texts for detailed information and recommended procedures.

Packaging

Product Name : Selenite Broth Base w/o Selenite

Product Code : DM402

Available Pack sizes : 100gm / 500gm

References

1. Klett A., 1900, Zeitsch Fer Hyg. Und. Infekt., 33: 137.
2. Guth F., 1916, Zbl. Bakt. I. Orig., 77:487.
3. Leifson E., 1936, Am. J. Hyg., 24(2): 423.
4. Murray P. R., Baron E. J., Jorgensen J. H., Pfaller M. A., Tenover F. C., Tenover R. H., (Eds.), 8th Ed., 2003, Manual of Clinical Microbiology, ASM, Washington, D.C.
5. Hartman P. A. and S. A., Munich, 1981, J. Food Pract., 44: 385-386.
6. Chattopadhyay W. and Pilford J. N., 1976, Med. Lab. Sci., 33:191.

Further Information

For further information please contact your local MICROMASTER Representative.



MICROMASTER LABORATORIES PRIVATE LIMITED




DM402PSS,QAD/FR/024,Rev.00/01.01.2018

Unit 38/39, Kalpataru Industrial Estate,
Off G.B. Road, Near 'R-Mall' , Thane (W) - 400607. M.S. INDIA.
Ph: +91-22-25895505, 4760, 4681. Cell: 9320126789.
Email: micromaster@micromasterlab.com



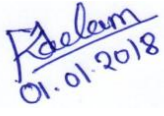

PRODUCT SPECIFICATION SHEET

Prepared By	Checked By	Approved By
 01.01.2018	 01.01.2018	 01.01.2018
Microbiologist	Head Quality Control	Head Quality Assurance

Disclaimer :

All Products conform exclusively to the information contained in this and other related Micromaster Publications. Users must ensure that the product(s) is appropriate for their application, prior to use. The information published in this publication is based on research and development work carried out in our laboratory and is to the best of our knowledge true and accurate. Micromaster Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are intended for laboratory, diagnostic, research or further manufacturing use only and not for human or animal or therapeutic use, unless otherwise specified. Statements included herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.

