



## PRODUCT SPECIFICATION SHEET

---

### Rabbit Plasma Fibrinogen Trypsin Inhibitor Supplement (MS118)

#### Intended Use

This supplement is recommended for enumeration of coagulase positive *Staphylococi*.

#### Composition

Each vial contains the following:

*Ingredients	Concentration
Bovine fibrinogen	0.375g
Rabbit plasma	2.50ml
Trypsin inhibitor	2.50mg
Potassium tellurite	2.50mg

#### Directions

1. One vial sufficient for 100ml medium.
2. Rehydrate the contents of 1 vial by adding 10 ml of sterile distilled water.
3. Turn the vial upside-down slowly two to three times avoiding frothing. Alternately, the contents may be dissolved by gentle stirring for 1-2 minutes using a sterile glass rod or a stirrer.
4. Immediate dissolution is not obtained. For better dissolution, keep it standing for 1-2 hours.
5. Slight particulate matter may be observed which does not affect the performance parameter of the supplement.
6. Aseptically add the supplement in 90 ml of sterile, molten Baird Parker Agar Base (DM034) / Baird Parker Agar Base (FPT) (DM1031) cooled to 48°C.
7. Add the supplement slowly from the side of the flask with simultaneous gentle swirling of the flask to mix the supplement and the medium homogeneously.
8. Pour into sterile petri plates and use immediately.

#### Storage and Shelf Life

Store at 2 - 8°C. Use before expiry date on the label.

#### Note: \* Not for medicinal use

The temperature of the molten, sterile Baird Parker Agar Base should be strictly 48°C at the time of addition of the supplement. Any variation in this temperature may give erroneous false negative results

#### Further Information

For further information please contact your local MICROMASTER Representative.



#### MICROMASTER LABORATORIES PRIVATE LIMITED


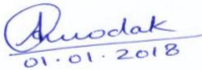

Unit 38/39, Kalpataru Industrial Estate,  
Off G.B. Road, Near 'R-Mall', Thane (W) - 400607. M.S. INDIA.  
Ph: +91-22-25895505, 4760, 4681. Cell: 9320126789.  
Email: [micromaster@micromasterlab.com](mailto:micromaster@micromasterlab.com)

MS118PSS,QAD/FR/024,Rev.00/01.01.2018.





## PRODUCT SPECIFICATION SHEET

Prepared By	Checked By	Approved By
 01.01.2018	 01.01.2018	 01.01.2018
Microbiologist	Head Quality Control	Head Quality Assurance

**Disclaimer :**

All Products conform exclusively to the information contained in this and other related Micromaster Publications. Users must ensure that the product(s) is appropriate for their application, prior to use. The information published in this publication is based on research and development work carried out in our laboratory and is to the best of our knowledge true and accurate. Micromaster Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are intended for laboratory, diagnostic, research or further manufacturing use only and not for human or animal or therapeutic use, unless otherwise specified. Statements included herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.

